

Caya® contoured diaphragm

intravaginal barrier device for her



PRODUCT BENEFITS:

- easy to insert and remove
- one size fits nearly every woman
- comfortable and practical for both partners

PRODUCT SPECIFICATIONS:

- elastic body: 50 shore-A durometer medical-grade silicone
- core: contoured, one-piece polymer spring
- overall length: 75 mm
- overall width: 67 mm
- membrane: 0.010" nominal thickness

PRODUCT FEATURES:

dynamic contours

- folds in optimal insertion shape
- embeds gently in supportive tissue

gentle dynamics

- soft and easy folding during insertion and removal
- slides gently out when removed

cervical cup membrane

- sized to fit a range of cervical positions
- surrounds the cervix

relief arch

- eases insertion or rather correction of the position
- avoids exertion of pressure on urethra

firm insertion edge


- folds compactly for easy insertion
- stable shape during insertion
- supports cervical cup in posterior fornix

grip dimples (nubs)

- enhanced grip when slippery
- cue for bending locations

removal dome

- enables easy hooking for removal (from top to bottom)
- enhanced insertion and positioning in anterior pubic bone

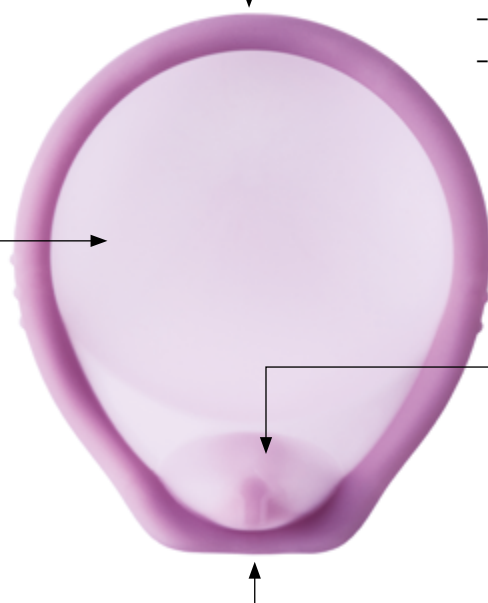

Each Caya® contoured diaphragm was separately checked for potential errors (100% control).

Produced and controlled in Germany. Further information at: www.caya.eu/en or www.caya.us.com



KESSEL MEDintim GmbH
Nordendstr. 82 - 84
64546 Mörfelden-Walldorf, Germany

Tel. (06105) 20 37 20
Fax (06105) 20 37 221
E-Mail: service@medintim.de

CE 0483

Medical device
Made in Germany



MEDintim
personal healthcare